
waterpik®

Waterpik® Sensonic® Professional Toothbrush

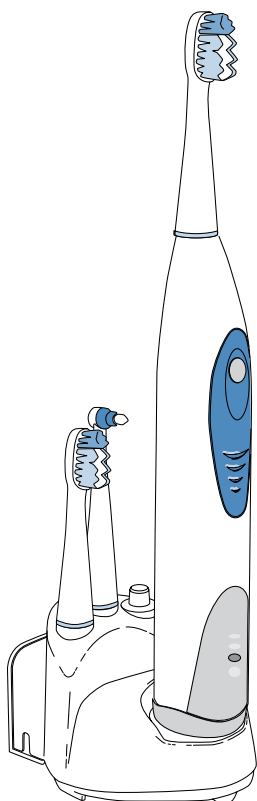**Waterpik® Sensonic®
Professional Zahnbürste**

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.....
www.waterpik.com

WELCOME

Congratulations, as an owner of the Waterpik® Sensonic® Professional Toothbrush, you are now on your way to better plaque removal, naturally whiter teeth, and a healthier smile!

The Sensonic® Professional Toothbrush utilizes an ultra high bristle speed and wide range of brushing motion. As a result, the Sensonic® Professional Toothbrush generates superior cleaning action and delivers improved plaque removal, healthier gums, and a fresh clean feeling.

Thank you for choosing the Waterpik® Sensonic® Professional Toothbrush.

Important Safeguards	4
Product Description	6
Getting Started	8
Care of your Waterpik® Sensonic® Professional Toothbrush	10
Limited Two-Year Warranty	11

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions should always be followed.

PLEASE READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle the plug with wet hands.
- Do not immerse charger in water or other liquid.
- Do not handle charger base while bathing.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not reach for a charger that has fallen into any liquid. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

WARNING:


To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has dropped into any liquid. Contact Waterpik International, Inc. or your local distributor which can be found at www.waterpik.com/intl.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Use only those attachments, including toothbrushes and interdental brushes, which are recommended by Water Pik, Inc.
- Keep charging unit and cord away from heated surfaces.
- Do not operate or use where oxygen or aerosol sprays are being used.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult with your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.


- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Be sure only the bristles of the toothbrush touch your teeth and gums when using the device.
- To ensure that your charger base does not accidentally fall and submerge in water, mount it on the wall.
- Pacemaker wearers and persons with medical implants containing magnets or electronics: contact your physician before using this product.

SAVE THESE INSTRUCTIONS

PRODUCT DESCRIPTION

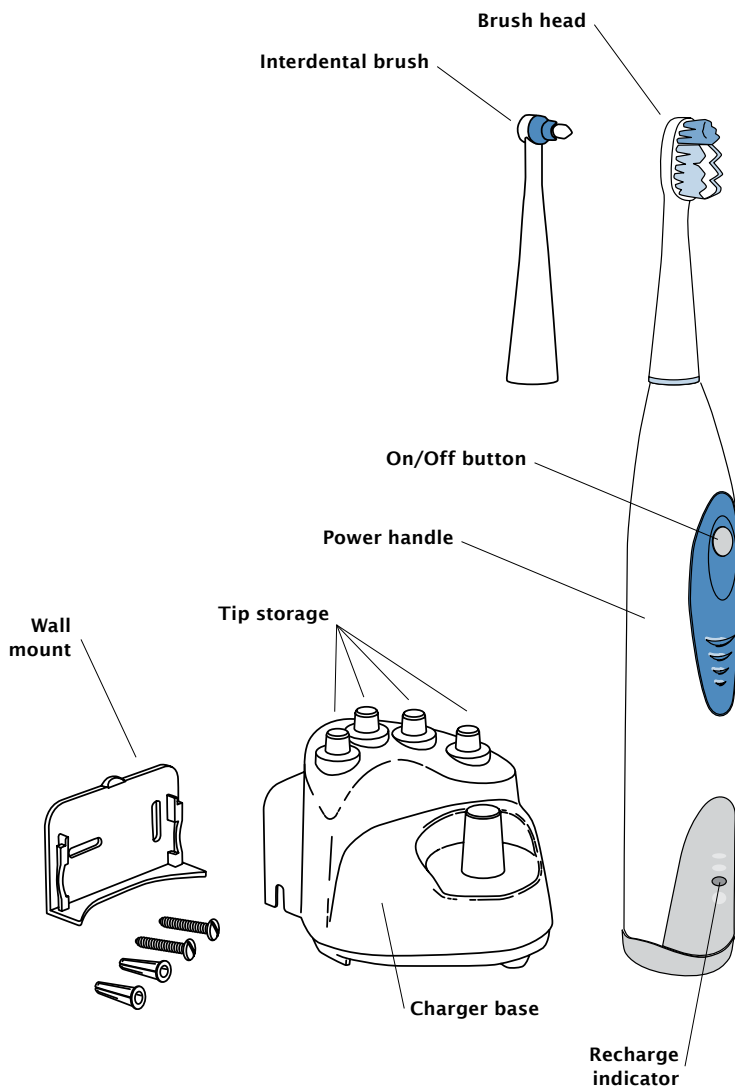
**Your Sensonic® Professional
Toothbrush can be easily mounted
to the wall using the enclosed
items.**

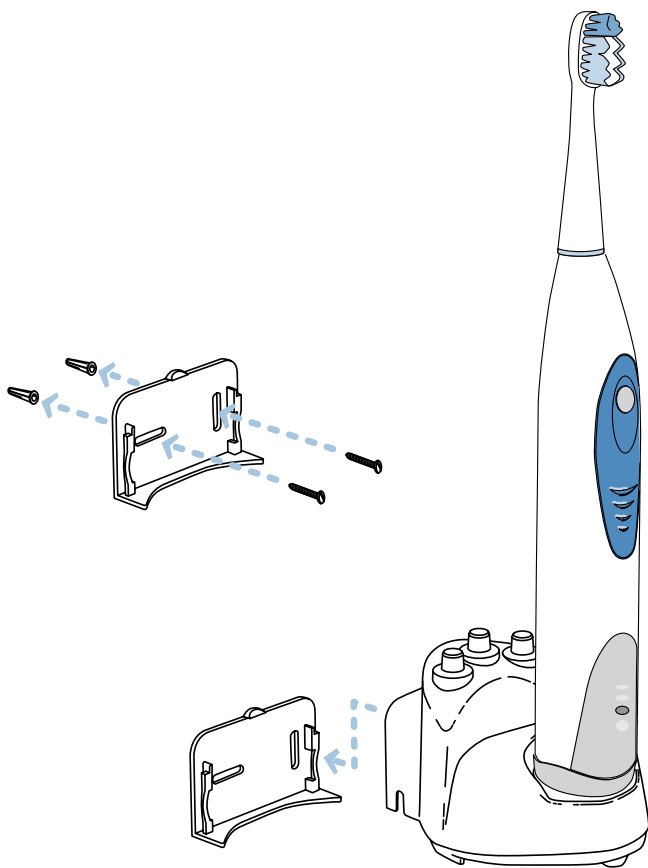GETTING STARTED


Charging the Sensonic® Professional Toothbrush **Charge for 16 hours before first use.**

The charge indicator on the handle will flash when charging and will stay on continuously when the unit is fully charged. A full charge will last for up to 1 week of typical use (2 minutes of brushing twice a day).

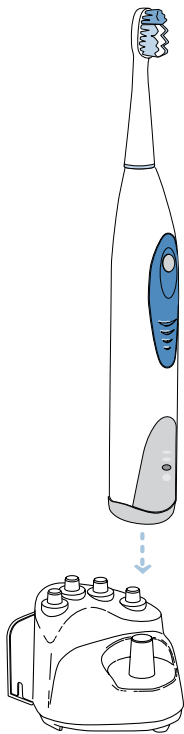
Sonic brushing

- Slide the brush head onto the handle. Push brush head down until it locks into place on handle.
- Place the brush head into your mouth before turning the unit on.
- Select low or high brushing speed by pressing power button after unit is on. Press button once for low speed, and twice for high speed.
 - Low speed – gentle cleaning and gum massage
 - High speed – superior plaque removal
- For best results, position brush along gum line at a 45-degree angle.
- Use a lighter brushing pressure than when using a manual toothbrush to optimize sonic cleaning effectiveness.

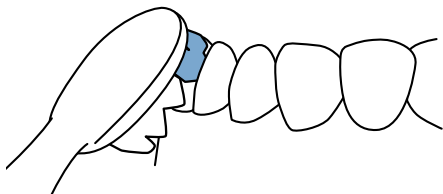- Move the brush slowly around entire mouth. Many dentists recommend spending 30 seconds in each of the four quadrants of your mouth.
- Quad Timer: your Sensonic® Professional Toothbrush will create a short stuttering sound every 30 seconds, to remind you to move to the next section of your mouth.

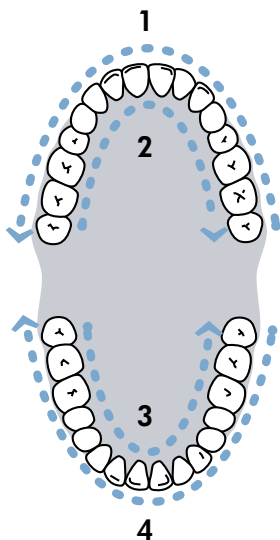Using the interdental brush

- Slide the interdental brush head onto the handle. Push brush head down until it locks into place.
- Brush may be used with or without toothpaste.
- Place brush tip perpendicular to gum line between two teeth.
- Turn on power handle and lightly follow crevices between teeth.

CARE OF YOUR WATERPIK® SENSONIC® PROFESSIONAL TOOTHBRUSH

1. After each use, turn off unit, remove the brush head from the handle and rinse brush head under water.
2. Wipe the handle dry after each use.
3. Periodically, clean the charger base by unplugging the unit and wiping with a damp cloth.

Service

There is no required routine service maintenance and this product cannot be repaired.

For all your warranty and accessory needs, please visit www.waterpik.com/intl. Please do not ship product to Waterpik International, Inc. offices, this will delay service. Refer to serial number and model number in all correspondence. These numbers are located on the bottom of your unit.

To locate your local distributor, go to www.waterpik.com/intl.

Battery Disposal

The unit contains a Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable. Batteries are permanently installed for performance reasons and are not removable for service or replacement. Batteries may be removed for recycling by dismantling the unit.

Important Precautions

Using Arm & Hammer® PeroxiCare® toothpaste, abrasive or chemical cleaners may damage the plastic housing of this product. Use of these products voids this product warranty.

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as brush heads.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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Arm & Hammer® PeroxiCare® are not registered trademarks of Church & Dwight Co, Inc.

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Canada, Chile, China, Czech Republic, Hungary, India, Mexico, New Zealand, Pakistan, Poland, Russian Federation, Singapore, South Africa, Switzerland, Taiwan, Ukraine, United Kingdom and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Canada, EU, Japan, Mexico, and the United States.

Sonicare® is a trademark of Water Pik, Inc. registered in Canada, EU, Switzerland, and the United States.

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1730 East Prospect Road
Fort Collins, CO 80553-0001
USA
1-800-525-2774

Europe Office:

Castle Court
41 London Road
Reigate, Surrey RH2 9RJ
United Kingdom

www.waterpik.com

Form No. 20012546-F AA

©2009 Water Pik, Inc.

644 Eden Park Drive
Cincinnati, Ohio 45202
t 513.961.6225 | f 513.281.3611

JOB INFORMATION

Product Description: WaterPik

Multilanguage Instruction Manual Model SR-1000

Client: WaterPik

Manufacturer:

Blueprint No: 1x67 164, 5x7 164d

UPC No:

Last Revision Date: 1/13/10

STUDIO INFORMATION

Job No: 2272.0006.028

File No: 2272.0006.028 SR 1000 Sansonic ML Guid

Studio Contact: Client Service: Mark Morse

Creative Lead: Scott Smith

Production Mgr: Doug Bennett

CLIENTS PACKAGING APPROVALS

BRAND	DATE
LEGAL	

IMPORTANT NOTICE

I/We have checked this artwork for accuracy. Final approval is the client's responsibility. Please double check for errors before reproduction. We recommend that color proofs be submitted to us for approval before production.

TECHNICAL INFORMATION

Application/Version: Adobe InDesign CS3

Layers/Pages: 5/1

Photos/Images:

Fonts Used: Refer to file

Separation: Name

Release Date:

Approved For Release _____ / /

COMPUTER EQUIVALENTS

TINTS / BUILT COLORS

brandimage APPROVALS

gmb	_____	_____	_____
Art Director	Date	Proofreader	Date
Account Executive	_____	_____	_____
	Date	Proofreader	Date
Design Director	_____	_____	_____
	Date	Proofreader	Date