waterpik®

Waterpik® Water Flosser plus Sonic Toothbrush Model WP-950/952

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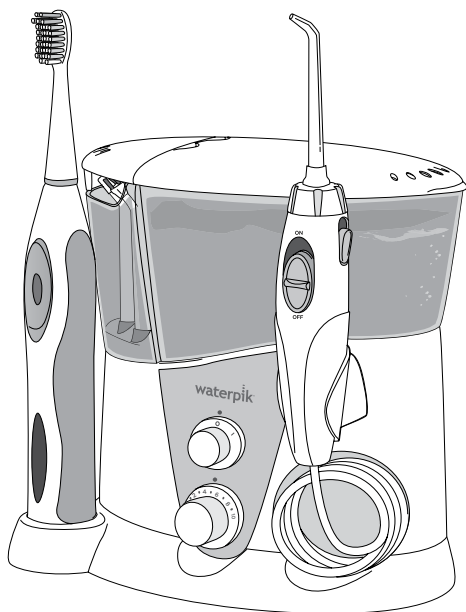
IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- If product is opened/disassembled for ANY reason warranty is **VOID**.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Water Pik, Inc. via email at InternationalService@waterpik.com
- Do not direct water from the Water Flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not use Water Flosser or Sonic Toothbrush outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

- Do not use iodine, or water insoluble concentrated essential oils in the Water Flosser. Use of these can reduce product performance and will shorten the life of the product.
- Remove any oral jewelry prior to use of the Water Flosser or Sonic Toothbrush.
- Do not use either the Water Flosser or Sonic Toothbrush if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist and/or medical professional before using the Water Flosser or Sonic Toothbrush or any other oral hygiene aid.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device containing magnets or electronics, contact your physician or the device manufacturer prior to using the Sonic Toothbrush.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- Do not use the Water Flosser for more than 5 minutes in each two-hour period.
- This product does not contain any user serviceable parts in the toothbrush handle or Water Flosser unit. Do not attempt to open either component except for battery disposal. See Battery Removal Section.
- Using toothpastes containing baking soda or peroxide, abrasive or chemical cleaners may damage the plastic housing of this product.
- Closely instruct and supervise children 8 years and above and individuals with special needs in the proper use of the Water Flosser and Sonic Toothbrush.
- Children should be supervised to ensure that they do not play with the appliance.
- There are no consumer serviceable parts in this appliance and it does not require maintenance.
- Fill Water Flosser reservoir with warm water or other dental professional recommended solutions only.


SAVE THESE INSTRUCTIONS.

TABLE OF CONTENTS

Product Description	5
Getting Started With Your Water Flosser	6
Cleaning and Troubleshooting Your Waterpik® Water Flosser	8
Getting Started With Your Sensonic® Professional Plus Toothbrush	9
Care of Your Sensonic® Professional Plus Toothbrush	10
Water Flosser and Sensonic® Professional Plus Toothbrush Limited Two-Year Warranty	11

PRODUCT DESCRIPTION


1. Brush head
2. Hinged tip/brush head storage lid
3. Tip/brush head storage
4. Toothbrush On/Off button
5. Toothbrush power handle
6. Toothbrush Recharge Indicator
7. Power cord
8. Power cord storage (under unit)
9. Toothbrush charger
10. Tip
11. Rotating direction control knob
12. Hinged reservoir lid
13. Water Flosser tip eject button
14. Reservoir
15. Water On/Off Switch
16. Water Flosser handle
17. Water Flosser power on/off switch
18. Water Flosser pressure control dial
19. Base

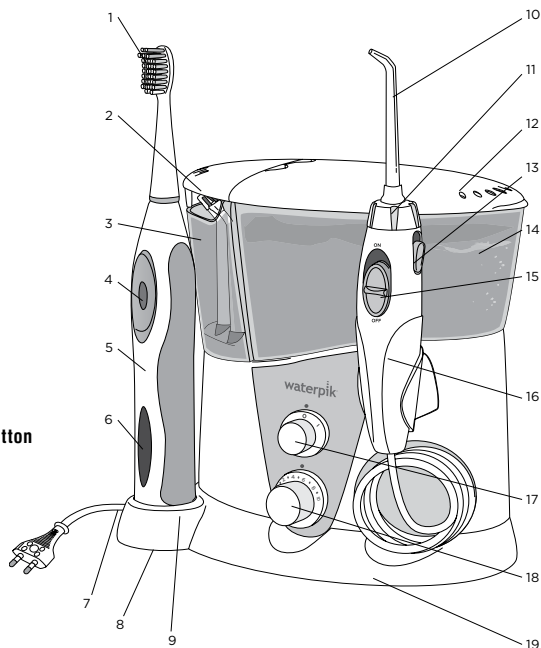TIPS

- A. Orthodontic Tip***
- Braces/General Use
- B. Pik Pocket™ Tip***
- Periodontal Pockets/Furcations
- C. Plaque Seeker® Tip***
- Implants/Crowns/Bridges/Retainers/General Use
- D. SensoNik® Large Brush Head***
- Plaque Removal
- E. SensoNik® Small Brush Head***
- Plaque Removal
- F. Classic Jet Tip****
- General Use
- G. Tongue Cleaner****
- Fresher Breath


NOT ALL TIPS INCLUDED IN ALL MODELS

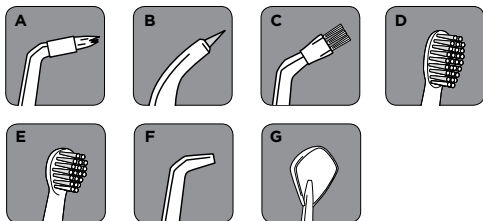*Replace Every 3 months

**Replace Every 6 months


GETTING STARTED

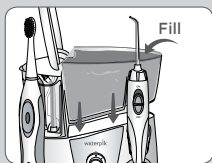**Powering the Unit**

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

Preparing the Reservoir


Lift the reservoir lid, remove the reservoir, and fill with lukewarm water. Replace the reservoir over the base and press down firmly and close the lid.

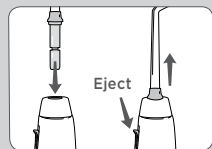**Inserting and Removing Tips**


Insert tip into the center of the Water Flosser handle, pressing down firmly until tip clicks into place. The colored ring will be flush with the top of the handle if the tip is correctly locked into place.

To remove tip from the handle push the tip eject button on the handle.

Caution: Do not press tip eject button while unit is running.

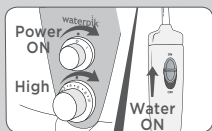**Prime Unit (First Time Use Only)**

Set pressure control dial on base to high (pressure setting No. 10), point the handle and tip into sink. Turn unit ON (I) and water ON until water flows from the tip. Turn Power Off (O).

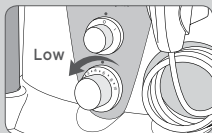**Adjusting the Pressure Setting**

Turn the pressure control dial on the base of the unit to the lowest setting (pressure setting No. 1) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

**Using the Water Flosser**

Always place the tip in your mouth before turning the water ON. Turn power on, lean over sink and slightly close lips enough to prevent splashing while still allowing water to flow from mouth into the sink. Turn water switch ON. When finished, turn water switch and power OFF. During use you can stop the flow of water by sliding the water switch on handle to the OFF position.

**Recommended Technique**

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

When Finished

Turn the unit OFF (0). Empty any liquid left in the reservoir. Unplug the unit.

NOTE: If you used mouthwash or other solution besides water, see "Using Mouthwash and Other Solutions" section.

NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

Orthodontic Tip




To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.


Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline. The toothbrush tip can be used with or without toothpaste. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes – much like you would with a manual toothbrush.

Plaque Seeker® Tip

To use the Plaque Seeker® tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

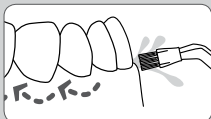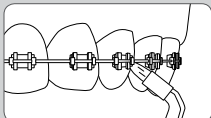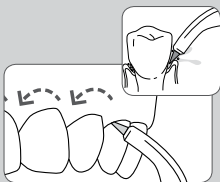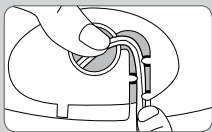Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Excess Power Cord Storage

After determining the length of power cord necessary, push extra cord into the cavity on the bottom of the unit and insert cord into the retention slots.

Cleaning

Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

The accessory storage container and the water reservoir are removable for easy cleaning and are top rack dishwasher safe.

Removing Hard Water Deposits

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 2 tablespoons of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every 1 to 3 months to ensure optimal performance.

Leakage

If leaking due to a damaged hose, please email InternationalService@Waterpik.com

Service Maintenance

Waterpik® Water Flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all your service needs email Water Pik, Inc. at InternationalService@Waterpik.com. Please do not ship product to Waterpik International, Inc. offices, this will delay service. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.


Still have questions? Email us at InternationalService@Waterpik.com.

Charging the Sonic® Professional Plus Toothbrush

Charge for 24 hours before first use. All 3 charge indicator lights will stay on continuously when the unit is fully charged. A full charge will last for approximately 1 week of typical use (2 minutes of brushing twice a day). If the toothbrush is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage.

Recharge Indicator Lights

Charge indicator lights will indicate level of battery charge for approximately 10 seconds after use and after removal from the charging base, and will then turn off.

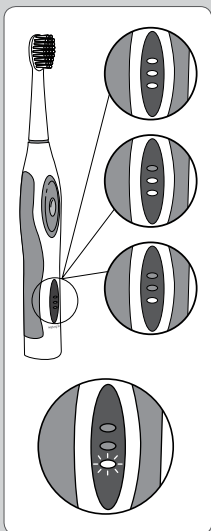Fully Charged

2/3 Charged

1/3 Charge life
(recommend recharge)



Unit will not operate if the bottom LED begins to blink rapidly.

In order for the toothbrush to begin to work again, it must be placed on the charger for approximately 1 hour.

CARE OF YOUR SENSONIC® PROFESSIONAL PLUS TOOTHBRUSH

Sonic brushing

- Slide the brush head onto handle and push it down firmly.
- Place the brush head into your mouth before turning the unit on.
- Press button once for low speed, and twice for high speed.
- If you hold the on/off button down briefly while in the low speed position the brush will shut off.
- For best results, position brush along gumline at a 45-degree angle.
- Use light pressure with a small, slow back and forth motion covering one or two teeth at a time.
- **Quad Timer:** Your Sensonic® Professional Plus Toothbrush will create a short pause every 30 seconds, to remind you to move to the next quadrant of your mouth.
- **2-Minute Timer:**
The brush will automatically shut off after 2 minutes of run time.

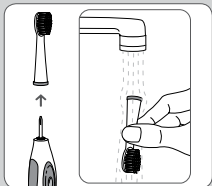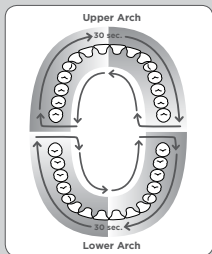Cleaning

1. **After each use,**
turn off unit, remove the brush head from the handle and rinse brush head under running water.
2. Wipe the handle dry after each use.

Service

Waterpik® toothbrushes have no consumer serviceable or replaceable electrical items and do not require routine service maintenance. For all your warranty and accessory needs, please email InternationalService@Waterpik.com. Please do not ship product to Waterpik International, Inc. offices, this will delay service. Refer to serial number and model number in all correspondence. These numbers are located on the bottom of your unit.

Disposal

This product contains rechargeable Nickel Metal Hydride batteries. The batteries are not replaceable.

Please do not dispose of the unit in the household waste. Remove the batteries per the instructions below and recycle at your local recycling center.

WATER FLOSSER AND SENSONIC® PROFESSIONAL PLUS TOOTHBRUSH LIMITED TWO-YEAR WARRANTY




Battery Removal

Note: This process is not reversible

Remove the battery only at the end of the product's life for recycling. Doing so earlier will destroy the product. Deplete the battery charge until the brush no longer operates. Observe basic safety precautions.

Remove the batteries as shown and recycle them and the rest of the product according to local regulations.

1. Use coin or washer to rotate cap counter-clockwise until cap protrudes.
2. Push down hard with shaft on hard surface.
3. Remove cap and pull internal body out. Cut tape and remove.
4. Lift battery pack up and cut battery wires, one at a time, to free from circuit board.

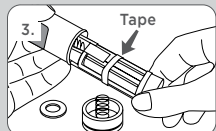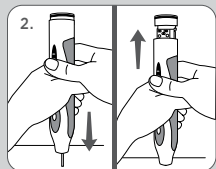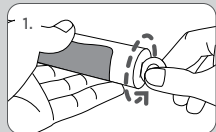Water Pik, Inc. warrants to the original purchaser/owner of this product that it is free from defects in materials and workmanship for two years from date of purchase. **Save your receipt as your proof of purchase date. Water Pik, Inc. reserves the right to request the product to be returned for analysis to confirm the claim falls within our warranty policy.**

We will replace any part of the product, which in our opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase and was used according to these instructions. For repair or replacement services provided for by the guarantee, please contact, via email, InternationalService@Waterpik.com. Upon request of the consumer, a copy of the guarantee will be provided by email or mail. This limited warranty excludes accessories or consumable parts such as tips, brush heads, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence.

Contact Water Pik, Inc. via email at InternationalService@Waterpik.com.

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waterpik®

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Brazil, Canada, Chile, China, Czech Republic, Denmark, EU, Finland, France, Germany, Hong Kong, Hungary, India, Israel, Italy, Japan, Kazakstan, Korea, Malaysia, Mexico, New Zealand, Norway, Pakistan, Philippines, Poland, Russian Federation, Saudi Arabia, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Turkey, Ukraine, United Kingdom, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Israël, Japan, Korea, Malaysia, Mexico, Philippines, Russian Federation, Saudi Arabia, South Africa, Switzerland, Taiwan, Trinidad & Tobago, Ukraine and the United States.

waterflosser® (stylized) is a trademark of Water Pik, Inc. registered in Brazil, Canada, China, Mexico, Switzerland and the United States.

Sensonic® is a trademark of Water Pik, Inc. registered in Canada, China, EU, Japan, Switzerland, and the United States.

Pik Pocket™ is a trademark of Water Pik, Inc.

Plaque Seeker® is a trademark of Water Pik, Inc. registered in Canada, China, EU, and the United States.

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Pik Pocket™ è un marchio di Water Pik, Inc.
Plaque Seeker® è un marchio di Water Pik, Inc. registrato in Canada, Cina, UE, e negli Stati Uniti.
Sensonic® es una marca comercial de Water Pik, Inc. registrada en Canadá, China, UE, Japón, Suiza y los Estados Unidos.
Pik Pocket™ es una marca comercial de Water Pik, Inc.
Pik Pocket™ is een handelsmerk van Water Pik, Inc.
Pik Pocket™ é uma marca comercial da Water Pik, Inc.